WiTscale S220

Bluetooth 4.0 Body Scale
With Body Fat Measurement

USER MANUAL


3. The Scale is equipped with auto-power off function. Once you finished the measurement, it will turn off in about 15 seconds.

CHANGE DISPLAY UNIT

- To change the display unit on the scale, push the button on the bottom of the scale. The display units will switch between lb, kg and st.

lb: pound
kg: kilogram
st: stone

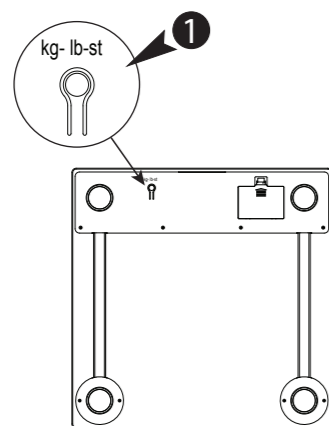Note:
Changing the display unit on scale will automatically change the measurement reading unit on Smartphone Apps. User can switch their desired display unit under setting tab on the Apps.

INTRODUCTION

Thank you for purchasing our product Bluetooth Body Scale S1(F) with Body Fat measurement. It supports iPhone 4S/ 5 and new generation iPad. We strongly recommend to read this Instruction Manual before you operate the product.

ABOUT WiTscale

WiTscale is a Bluetooth 4.0 connected digital body scale. WiTscale can seamlessly transmit your personalized weight, Body Fat and BMI data to your iPhone or iPad. All recorded data can be tracked and shared with your friends on facebook. Setting up your personalized profile and fitness target will take less than three clicks.

CAUTION

- Persons with pacemaker should not use this product
- Persons with metal implants will be unable to obtain accurate readings

BOX CONTENTS

- Your WiTscale will contain:
- WiTscale body scale unit
 - 3 AAA batteries
 - This user manual
 - Quick start guide

APPS


Setting up the body scale to work with your iPhone or iPad, user needs to download our free Apps from Apple AppStore at www.apple.com/apps and search keyword for "WiTscale".

Users can use the App to measure, track, record and share their weight and body fat information to extend the scale functions. User can still use the scale as stand alone device without the Apps but no body fat information display on the scale.


Note:
Currently, our WiTscale only support Bluetooth 4.0 iOS device. Older generations of Bluetooth technology is not supported.

WARNING INDICATION

- Weight exceed limit. User must step off immediately
- Low battery. User should replace the batteries for proper operation


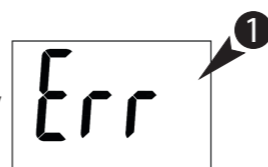
LCD DISPLAY

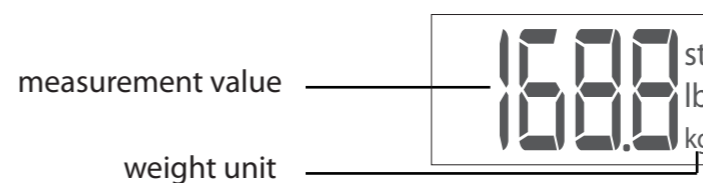BATTERY INSTALLATION

- Remove the battery door located on the bottom of the body scale. Insert 3 pieces of AAA batteries in the compartment with correct polarity indicated on the scale.
- After inserted batteries, the LCD will show full display for about 2 seconds
- Then the LCD will show 0.0 value and it is ready to weigh your body

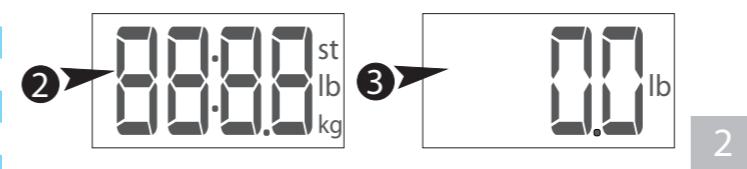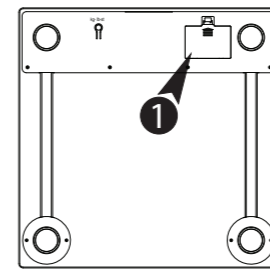BATTERIES

When the warning indication display LO, it is necessary to change the battery. On the bottom of the unit, open the battery compartment and replace with new battery only. Keep in mind, the display unit will be defaulted to LB (pound) after power reset.

ATTENTION


- Weight capacity: 7KG – 150 KG. Please weigh within this weight range for the use of the scale correctly and without causing damages.
- DO NOT place anything on the scale whilst not in use.
- This scale is to be used to measure the body weight of a human being. DO NOT use for other purposes.
- For cleaning of the scale, please use a soft fabric cloth with NO cleaning agents.
- Always stand in the centre of the scale for a correct reading and DO NOT operate scale whilst wet.
- For a precise result, only use scale on flat, hard & even surfaces.
- Store the scale in a cool dry place. DO NOT expose to extreme temperature or moisture.
- NEVER jump or bounce on the scale. DO NOT shake or try to disassemble the scale.
- Do not press or push hard on LCD display.

WEIGHT AND BODY FAT MEASUREMENT

Step on the Scale with bare foot and keep contact with electrodes (stainless surfaces). The weight will be displayed on LCD. Once the reading is steady, the display will flash 3 times to indicate the measurement is completed.

To make a new measurement, just come down from the scale and step on it again.

Put your feet properly on the scale as follow:

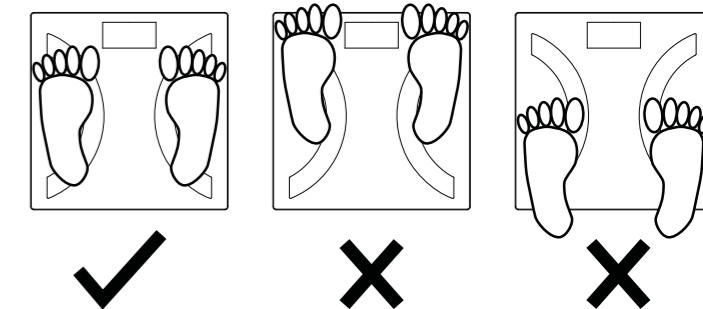- Notes:
- To obtain accurate measurement, make sure your feet can be closely contacted on all four stainless steel plates and stand straight/upward of your body.
 - The body scale should not be placed on any uneven or carpeted flooring to avoid unbalanced body.
 - To ensure accuracy, we suggest to use the scale at the same time every night (19:00 - 21:00)

SPECIFICATIONS

Weight Limits	: 5 to 150 Kg
Body Fat Limits	: 5 - 50 %
Power Supply	: AAA battery x 3
Operating Temperature	: 5 to 35°C
Wireless Standard	: Bluetooth 4.0
Frequency Range	: 2400 – 2480 MHz
Operation Range	: 10 meters

LIMITED WARRANTY

If WiTscale Digital Scale does not work properly because of defects in materials or workmanship within 90 days from the date of original purchase and when used in according to user manual provided with the scale. The option could either (a) repair your product with new or refurbished parts, or (b) replace with a new or refurbished unit. The final decision to repair or replace will be made by the US based distributor. All return or replacement postage will be paid by you. An original purchase receipt is required to obtain warranty service.
Important: This limited warranty DOES NOT COVER normal wear, result from accidents, damage in transit, misuse, abuse, mishandling and damage in humidity or liquids. Opening the hardware unit will automatically void the limited warranty.

FCC ID: PXM-S1

- This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and
(2) This device must accept any interference received, including interference that may cause undesired operation.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Any unused scale and batteries must be disposed of according to the local EPA regulations at the end of their service life.

Distributed by:
Season International Corp.
Park Ridge, Illinois, 60068
www.seasonations.com
Made in China

email support:
support@seasonations.com
toll-free support no.:
1-(866) 899-8924

iPhone and iPad are the registered trademarks of Apple Inc.

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iPhone e iPad son marcas registradas de Apple Inc.